


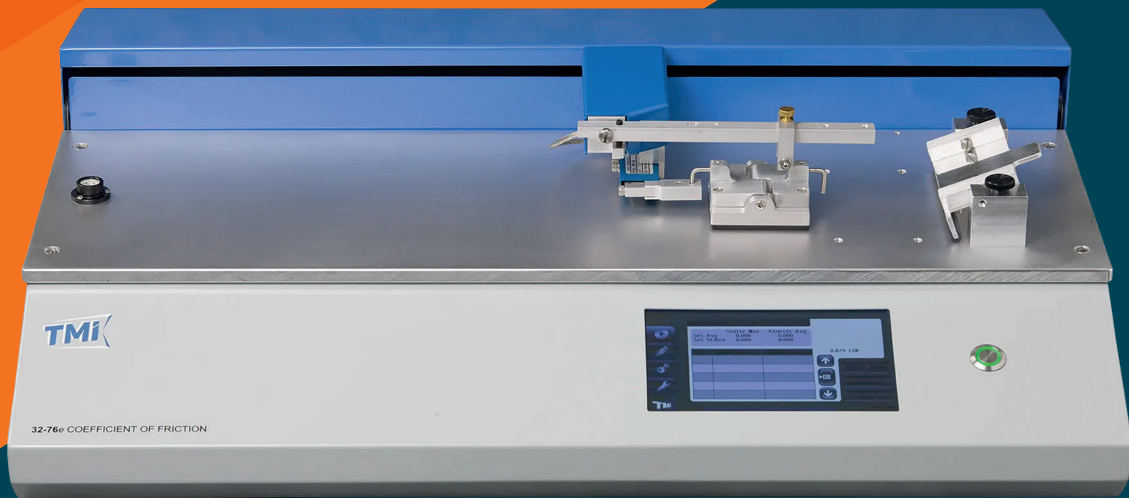

Coefficient of Friction/ Peel Tester

Model 32-76e

TMI Büchel is globally recognized as the leading supplier of coefficient of friction instrumentation.

Our NEW 32-76e coefficient of friction/peel tester uses advanced digital force signaling and high-speed data acquisition software to provide unmatched precision and repeatability in COF and peel testing.

Improved features include a color touchscreen display and intuitive software user interface for easy navigation and test method storage.

New peel testing capabilities include 180° peel and T-peel measurements.


Seal Strength/Peel Software

A comprehensive peel strength and heat seal application is included with the instrument. The operator can enter a variety of peel strength parameters designed to support most international peel test and heat seal international standards. Selectable units, pre-peel and peel distance settings allow flexible test method set-up and storage.

Results include average, maximum, minimum and standard deviation of an individual peel test. The 32-76e also provides result and test curve data output to GraphMaster Plus, PC-based software for data collection, test curve analysis and results storage. A selection of test fixtures is available to meet most peel applications.

T-Peel Fixture

is designed to measure the adhesive force of a heat seal or adhesive bonded substrate. The

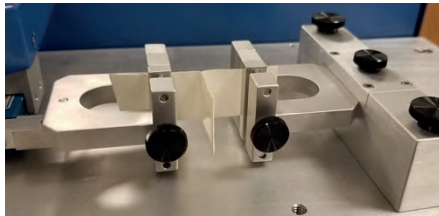

T-Peel Fixture (32-76-03)


fixture attaches directly to the moving arm of the 32-76 Friction/ Peel tester. A second clamp attaches directly to the instrument base. Included is one set (2) of test clamps, 30.5 mm (1.5 inch) wide, capacity is 100N.

Meets the following standards:

- ASTM D1876 T-Peel Test
- ASTM F 88 Seal Strength of Flexible Materials

180-Degree Peel Fixture

measures the release strength of a label or peeling strength of an adhesive

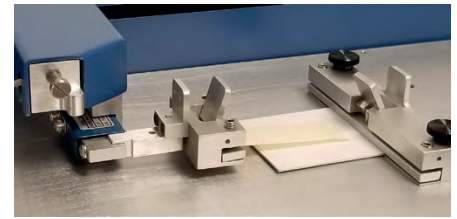

180-Degree Peel Fixture (32-76-04)


substrate. The fixture includes one set (2) of test clamps. One clamp is attached to the load cell with a clamping width of 50.8mm (2.0 inch) wide, maximum capacity 100N. A second clamp is attached to the instrument base, maximum sample width is 125 mm (5 inch). With this fixture, release strength of labels and tapes can be performed. Typically, the label or tape is attached to a rigid backing.

Meets the following standards:

- PSTC-101 International Standard for Peel Adhesion of Pressure Sensitive Tape
- ASTM D 3330 Standard Test Method for Peel Adhesion of Pressure-Sensitive Tape
- FINAT1- Peel adhesion 180°

Heated Platen Option


During the manufacture and converting process of plastic films, the material is subjected to elevated temperatures. When using a heated platen, it is common that the heated surface for most

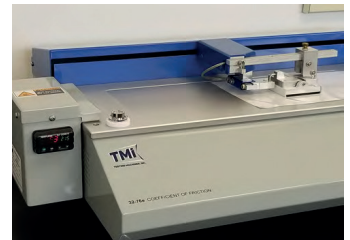

Heated Platen (32-76-08)

plastic films will produce an increase in COF. Normally a temperature increase will influence the friction of the packaging material. However, not all films will react in the same way to an increase in temperature. To measure this behavior, the 32-76e Friction/Peel tester is available with an optional heated platen. The temperature can be set from ambient to 162° C (325° F). The heated platen option can easily be installed in the field on an existing 32-76e unit (Part # 32-76-08). The heated platen option meets ASTM D 1894 standards for friction testing at heated temperatures.

GraphMasterPro™

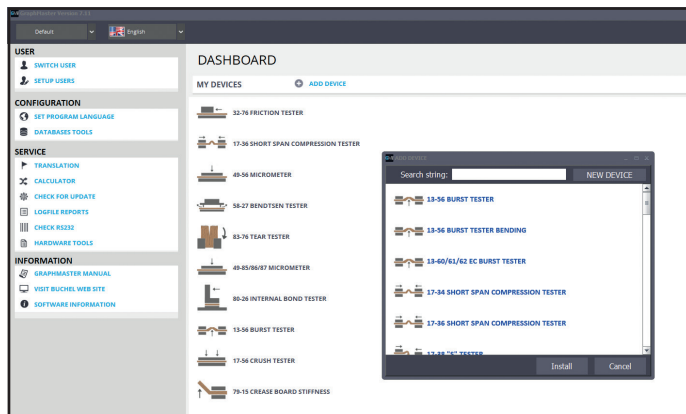


GraphMasterPro™ is a Windows™ based software data collection program capable of providing test curve analysis, result data, data storage, exporting and printed reports. The software automatically captures individual test results and/or raw data.

In addition to collecting data from instruments manufactured by TMI Büchel, GraphMasterPro™ can be configured to communicate with instruments manufactured by other instrument suppliers.


Part Number	Description
32-76-00-0002	32-76e Friction/Peel tester 110/220
32-76-00-0003	32-76e Friction/Peel tester w/Heated Platen 110/220
03-03-00-0001	GraphMasterPro™ software,
32-76-02	"B" Type Sled, 200 g, Magnetic Sled
32-06-02	"B" Type Sled, 200 g
32-76-03	T-Peel Accessory for 32-76e
32-76-04	180 Degree Peel Accessory for 32-76e
32-76-08	Heated Platen only
32-07-10	Calibration weight set 10-1000g
32-76-05	10N Load Cell
32-76-06	20N Load Cell
32-76-09	100N Load Cell
32-76-10	5N Load Cell

Test Method Storage

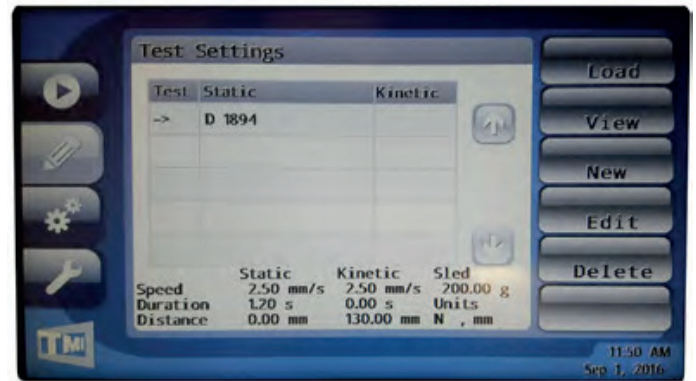

Test method storage allows users to save setup specifications, ensuring proper settings.

Single Screen Results

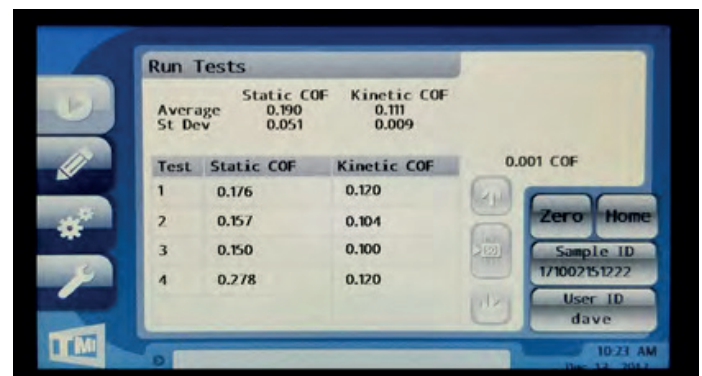

Test screen displays static, kinetic results and statistics on one screen, improving testing efficiency.



Magnetic "B" Type Sled, 200 g

Features

- 7" full-color digital touchscreen display
- High speed data collection and analysis for precise measurement of static COF with 500 readings during the first second
- Selection of test type-COF/Friction, T-Peel and 180° Peel
- Heated Platen Option, ambient to 162° C (325° F)
- Selectable units (COF, g, N, kg, lbs, and ounce)
- Selectable load cells from 5 to 100 N
- Automatically reports static and kinetic friction results after measurement
- Settable distance and time for static, kinetic and peel
- Peel results include average, minimum, maximum and SD of peel force
- Test storage including individual results up to 100 readings, average and standard deviation
- Direct drive arm with unique skid control.
- Sled-connecting mechanism ensures level pulling action
- New magnetic sled simplifies testing and throughput
- Direct output to PC-based GraphMasterPro™ software providing friction/peel curve analysis, storage and reports
- Mini USB output for PC connection (GraphMasterPro compatible)
- USB output for data storage

Applications

Plastics film, packaging, paper, labels, coatings, cartons, adhesives, foils and printed materials

Standards

COF conforms to ASTM D 1894, ISO 8295, TAPPI T 549

Peel conforms to ASTM D 1876, ASTM F 88, ASTM D 3330 PSTC-101, FINAT1

Technical Specification

Model	32-76-00-0002
Load cell	Interchangeable load cells: 5N, 10N, 20N, 50N and 100N
Speed Range	10-50 cm/min (4-20 in/per min)
Travel distance	2.5 to 30.5 cm (1 to 12 in)
Accuracy	+/- 0.1 % of load cell capacity
Resolution	0.001 gram based on 1000 gram/10N load cell (NOTE: Resolution will change depending on load cell in use)
Available Sleds	100g - 1816g (System software accepts variable sled weights) Please specify sled.
Weight	25kg (55lb)
Instrument Size	W x D x H 515mm x 495mm x 508mm (20.3in x 19.5in x 20in)
Power Input	120/240VAC, 50/60Hz
Data Output	USB flash drive, USB ESC/POS printer, mini USB, GraphMasterProTM

CoF Standards

TAPPI T549	Coefficients of Static and Kinetic Friction of Uncoated Writing and Printing Paper by use of the Horizontal Plane Method
TAPPI T816	Coefficient of Static Friction of Corrugated and Solid Fiberboard (Horizontal Plane Method)
ASTM D1894	Standard Test Method for Static and Kinetic Coefficients of Friction of Plastic Film and Sheeting

Peel Standards

ASTM D 1876	Standard Test Method for Peel Resistance of Adhesives (T-Peel)
ASTM D 3330	Standard Test Method for Peel Adhesion of Pressure-Sensitive (180° angle)
ASTM F 88	Standard Test Method for Seal Strength of Flexible Barrier Materials
PSTC-1016	Peel Adhesion of Pressure Sensitive Tape (180° angle)
FINAT 1	FTM 1 Peel Adhesion (180°) at 300mm per minute
EN 1939	Peel Adhesion Self Adhesive Tapes (180° angle)

Alternative Slide Angle Model

32-25-00	Incline Plane Coefficient of Friction Tester (Meets TAPPI T-815 and ASTM D202 with optional blocks)
-----------------	---

Contact Details

web. www.industrialphysics.com

email. info@industrialphysics.com

Disclaimer

The information contained in this document is liable to modification from time to time in the light of experience and our policy of continuous product development. Check the Industrial Physics website for the latest version.